

Aerogel RICH & TOP counter for super KEKB


Y.Mazuka Nagoya University

June 14-16, 2006
The 3rd SuperB workshop at SLAC

Introduction



Current BELLE performance

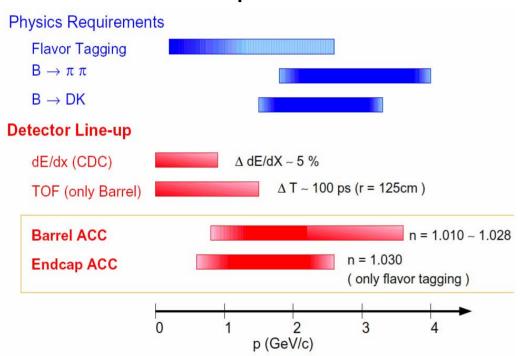

We cannot PID at high momentum region in the forward endcap

Present endcap-ACC is used only for flavor tagging



Further Improvement on π/K separation with the start of super KEKB


PID Target: π/K separation > 4σ @4GeV/c



Upgrading BELLE Detector



Two new particle ID devices, both RICHes

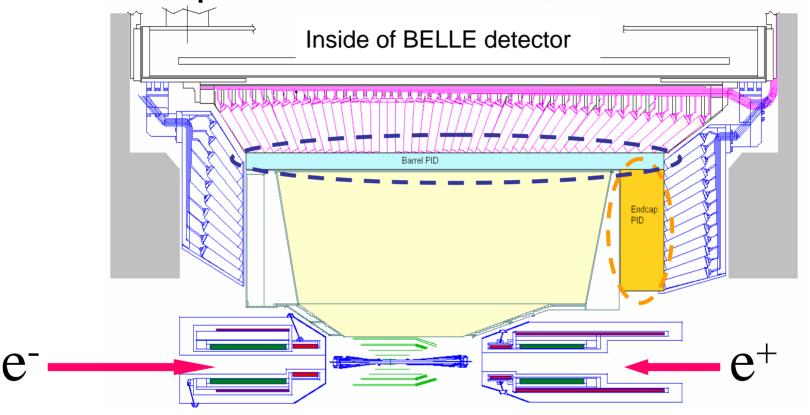


Endcap: Proximity Focusing Aerogel RICH(A-RICH)

Barrel: Time of Propagation Counter(TOP)




Outline



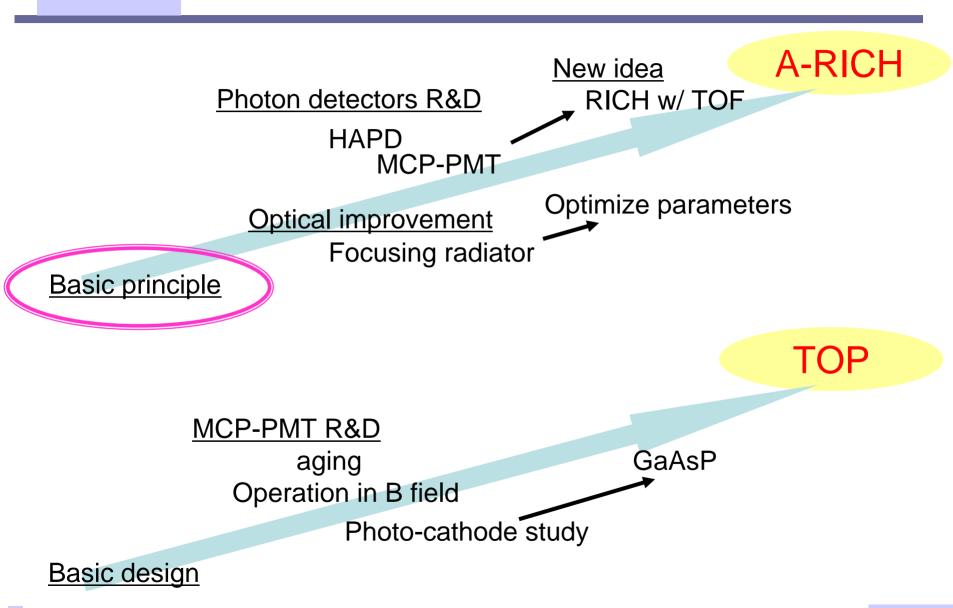


Endcap: Proximity Focusing A-RICH



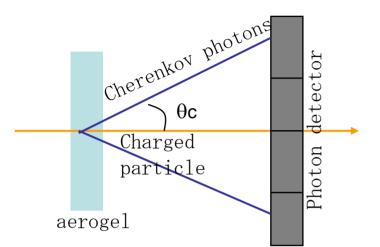

- PID by Cherenkov ring image emitted from aerogel radiator
- Cherenkov angle π/K (4GeV/c, n=1.05)

$$\theta_c^{\pi} - \theta_c^{K} \cong 23mrad$$

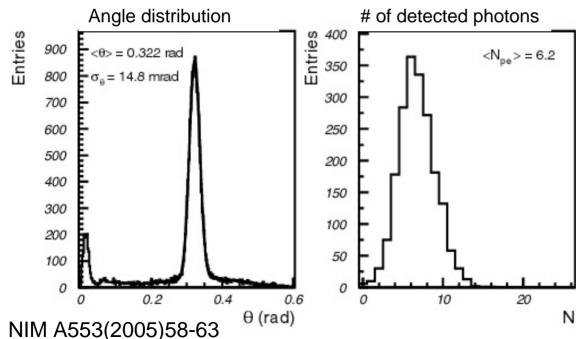

Typical beam test results π -, 4GeV/c, 2cm thick aerogel

 π/K separation ~ 4.0 σ It seems to reach our target, but


Can we improve the separation over 4σ ?

Focusing configuration




How to increase the number of photons without degrading the resolution?

Use radiator with gradually increasing refractive index in down stream direction

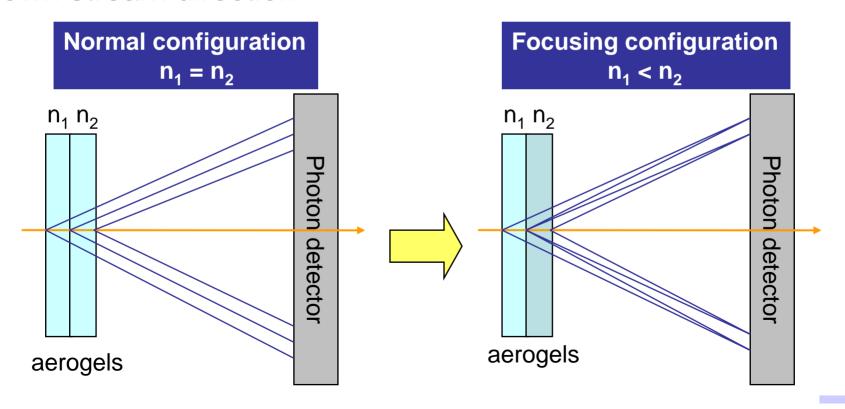

Results of focusing configuration



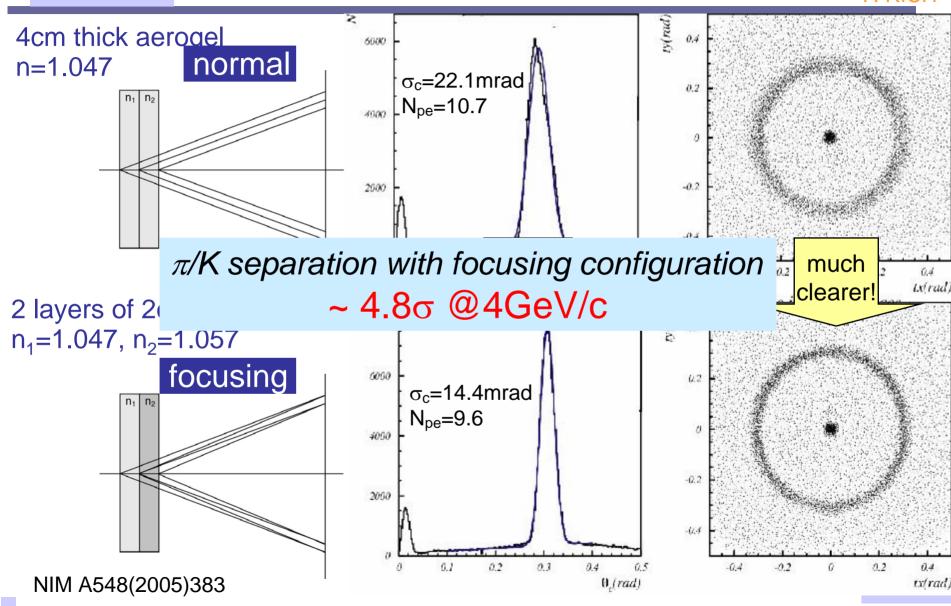


Optimization of dual radiator indices



Upstream aerogel: n₁=1.045


Downstream aerogel: n₂ is changed

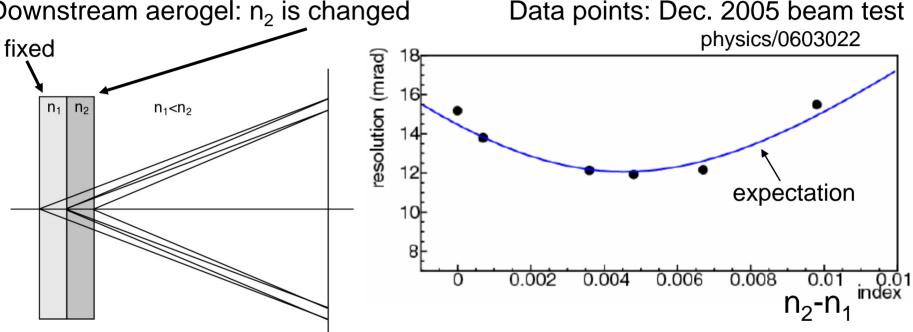

- Measured resolution is in good agreement with expectation
- Wide minimum region allows some tolerances(~0.003) in aerogel production

Outline



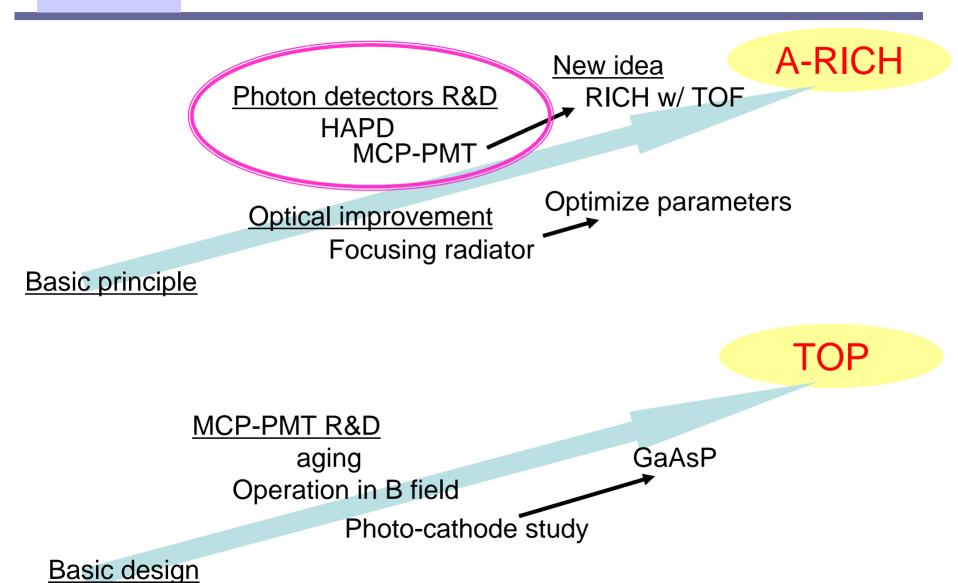

Photon Detectors for A-RICH



Requirements

- Working in B=1.5T
- Pixel size ~5-6mm
- Good sensitivity to single photon
- Large effective area
- Candidates
 - HAPD with large effective area

Y. Mazuka

MCP-PMT

Photon Detectors for A-RICH; HAPD




demerits

- Low gain (~10⁴)
- High noise rate

<u>merits</u>

- High efficiency
- High energy resolution

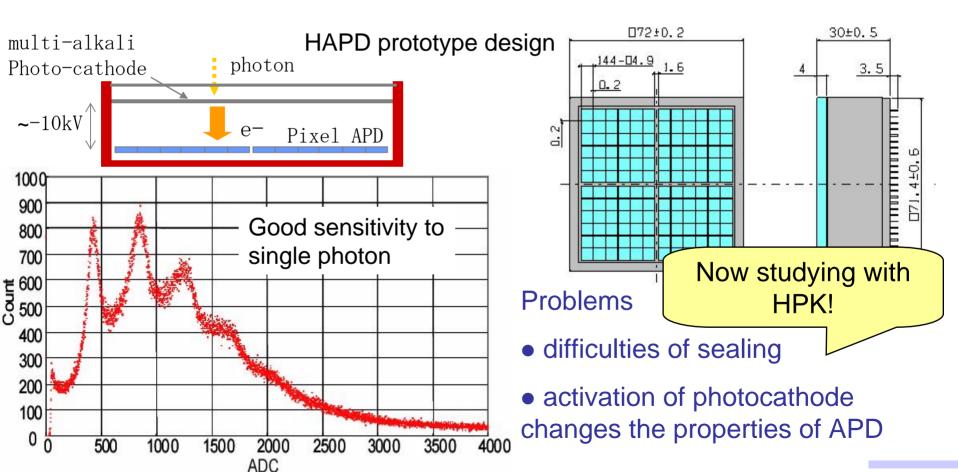


Photon Detectors for A-RICH; MCP-PMT



BURLE 85011 MCP-PMT

photo- cathode	Bi-alkali
MCP	25μm pores, 2 MCPs
gain	~ 0.6×10 ⁶
collection efficiency	~ 60%
dimensions	~ 71mm square
# of channels	8×8
pitch	~ 6.45mm
active area	~52%

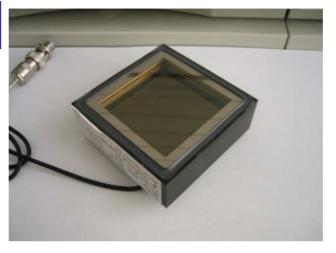

demerit

active area

merits

- High gain
- Good time resolution


TTS~50psec(single p.e.)
Can we use this merit?

A-RICH with TOF using MCP-PMT



New idea!

Fast photon detector enables A-RICH to have TOF info.

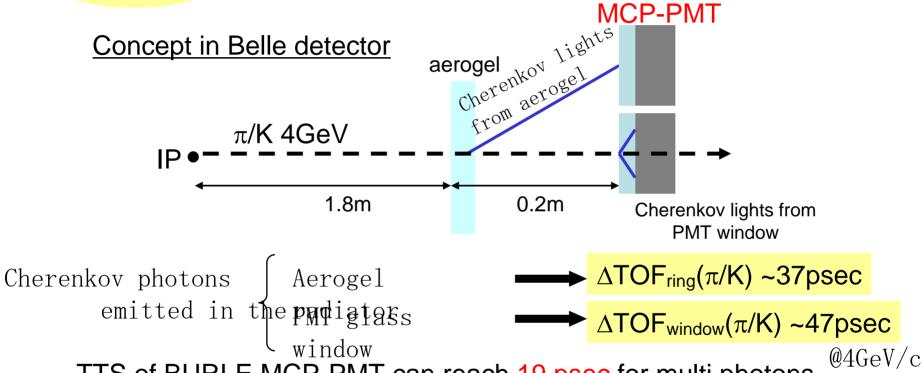

TTS of BURLE MCP-PMT can reach 19 psec for multi photons

 $S_{TOF}(\pi/K) > 2.4 \sigma$ for multi photons?


PID can extend lower momentum region

 $\Delta TOF(\pi/k)$ is bigger



A-RICH with TOF Beam test results



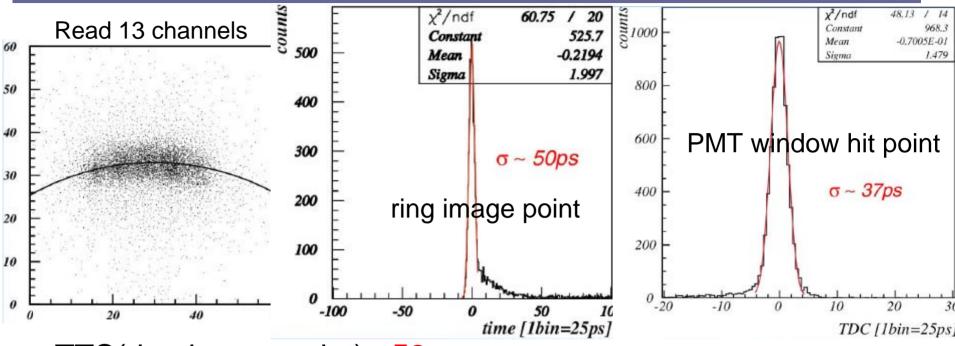

- TTS(ring image point) ~50psec
 - Time resolution per one track (N_{pe}~10) ~20psec
- TTS(PMT window hit point) ~37psec

 π/K separation with MCP-PMT $S_{TOF} \sim 2.2\sigma$ @4GeV/c




Consistent to

expectations

A-RICH with TOF PID at low momentum



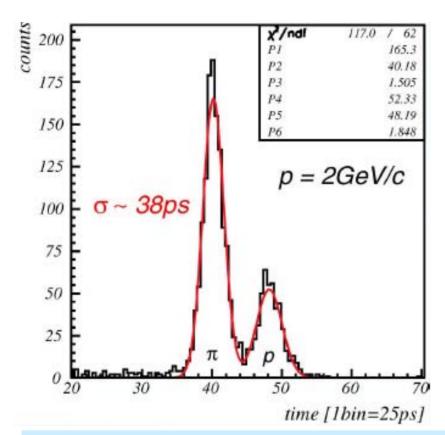

TOF test with pions and protons at 2GeV/c

Photons from PMT window π/p are well separated


Even in distance between start counter and MCP-PMT is 65cm, instead of 2.0m in Belle

At this test, π/p separation with MCP-PMT $S_{TOF} \sim 4.8\sigma$ @2GeV/c

A-RICH with TOF using MCP-PMT looks very promising

Outline




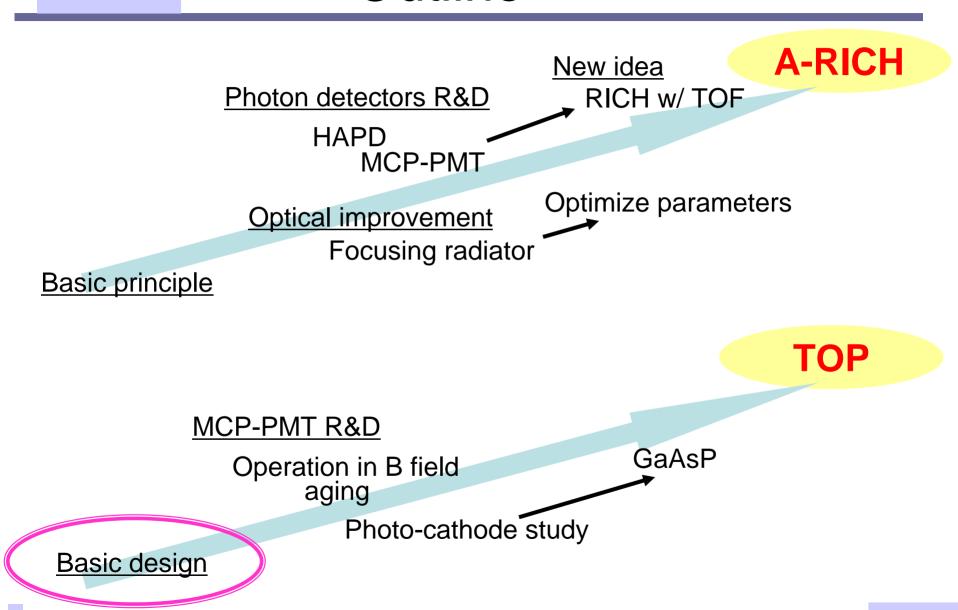


Barrel: TOP counter



Cherenkov ring imaging is used as timing information

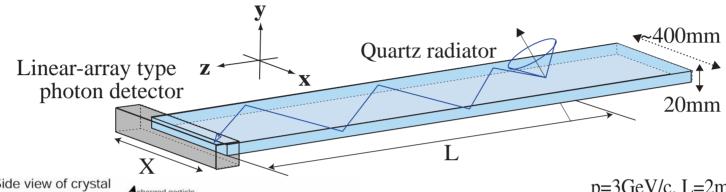


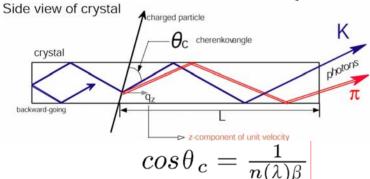

Difference of path length

→ Difference of time of propagation (TOP) (+ TOF from IP)

With precise time resolution (σ ~40ps)

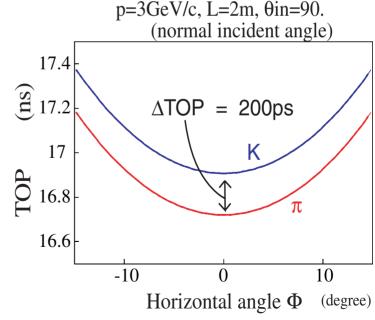

Photon Detector for TOP; MCP-PMT



- Requirements:
 - Good sensitivity to single photon
 - TTS~30ps (single photon)
 - working in 1.5T



- BURLE (25µm pores)
- BINP (6μm pores)
- HPK (6 and 10μm pores)




B=0T: all samples have good TTS(~30ps)

B=1.5T: BINP and HPK samples have

high gain($\sim 10^6$) and good TTS(~ 30 ps)



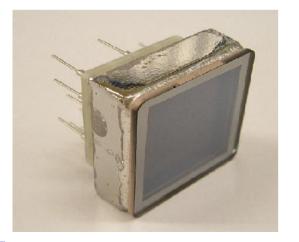


These samples were round shaped (1ch.)



MCP-PMT aging



Study of tubes w/ and w/o Al layer

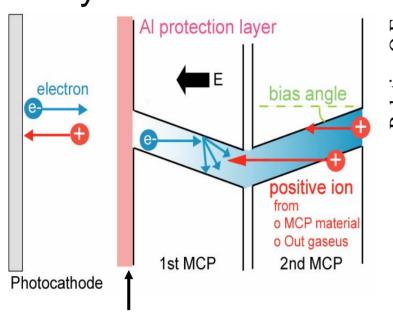

 $\lambda = 400 \text{nm}$ Relative Q.E Photo-cathode: multi-alkali 0.8 0.6 HPK w/ Al HPK w/o Al • 🔻 Russian w/ Al 0.4 o△ Russian w/o Al 0.2 20 Integrated irradiation (x10¹³photons/cm²) Time in Super B-factory (year)

Prevent feedback ions from reaching the photo-cathode

(It reduces collection efficiency by 60%)


HPK w/ Al survives over 13 years of operation!
Al layer is necessary


MCP-PMT with GaAsP



Expected performance

bialkali photo-cathode: π/K separation at 4GeV/c < 4 σ → chromatic dispersion

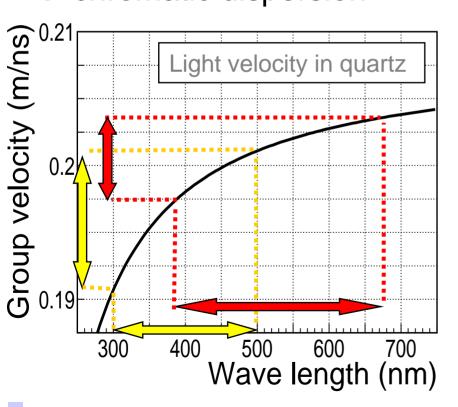

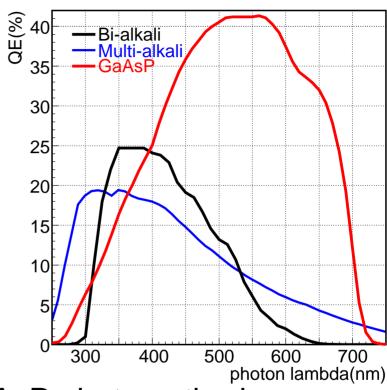

GaAsP photo-cathode:

- Higher Q.E.
- At longer wavelength
 - →less dispersion

 π/K separation > 4σ


The 3rd SuperB workshop, SLAC


GaAsP MCP-PMT development



- Square-shape MCP-PMT with GaAsP photocathode is under development with HPK
- First prototype
 - The same type as previous tubes
- Performance test
 - Gain
 - Time resolution



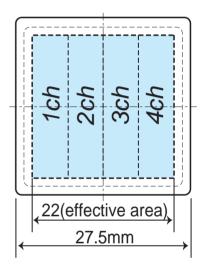


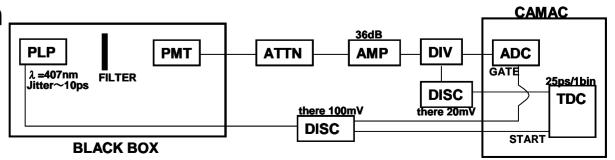



GaAsP MCP-PMT performance



Wave form, ADC and TDC distributions

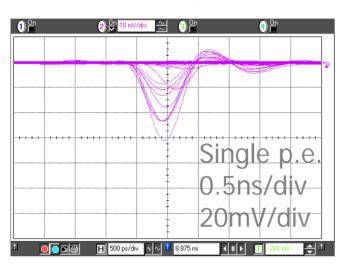

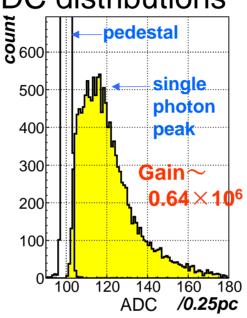

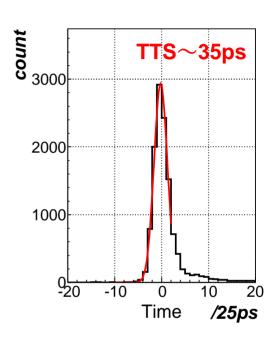

- Enough gain(~10⁶) to detect single p.e.
- Good time resolution (TTS~35ps) for single p.e.
- Next
 - Check the performance in detail
 - Life time of GaAsP photo-cathode tube

Summary



We are studying new types of RICH for super KEKB

Aerogel RICH counter for endcap

- Test the focusing configurations
 - We studied about optimal parameters
- More studies: RICH with TOF (using MCP-PMT)
 - Extend PID ability into low momentum region

TOP counter for barrel

- N Both RICHes(A-RICH, TOP) look very promising
 - π/K separation can be over 4σ @4GeV/c

The 3rd SuperB workshop, SLAC

- But there is still a lot of work to be done!
 - It will reduce the effect of chromatic dispersion




23/24

Tasks for practical use



A-RICH

- Photon detectors
 - Develop HAPD & MCP-PMT in parallel
- Readout system
 - ASIC
- Mechanical design
 - Line up of photon detectors and radiators

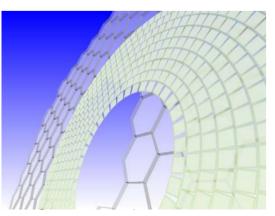





- Make practical tube
- Aging of tube with GaAsP
- Readout system
 - TAC
- Test of prototype
 - Line up of photon detectors and radiators




Backup




25/24

Optimal aerogel thickness




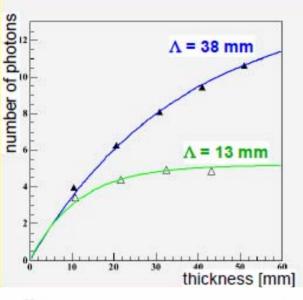

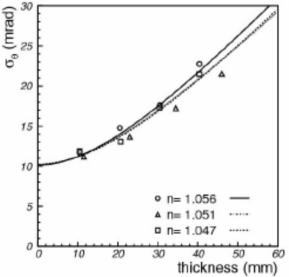

Cherenkov angle resolution per track

$$\sigma_{track} = \frac{\sigma_{\theta c}}{\sqrt{N_{pe}}}$$

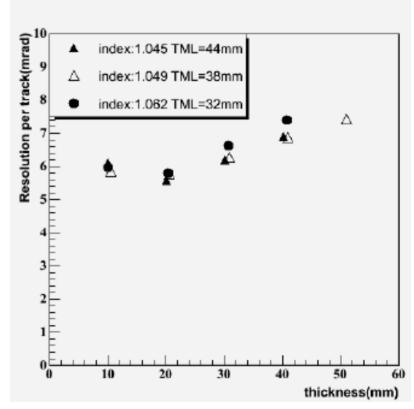


The best resolution ~ 5.5mrad at 2cm thick aerogel

RICH with TOF





 Δ TOP(K- π) \sim -5psec (Difference of TOP from emission point of Cherenkov light to PMT)

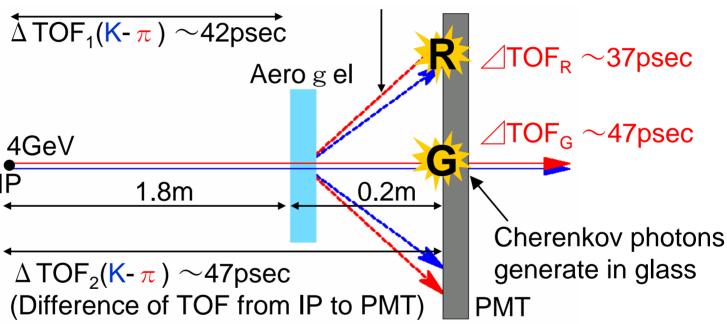

Time resolution of 10 psec has been achieved with HPK MCP-PMT @ Nagoya university.


Time resolution of BURLE MCP-PMT can reach 19 psec for multi photons.

More than 2.4σ for multi photons?

A-RICH with TOF using MCP-PMT



Setup of beam test

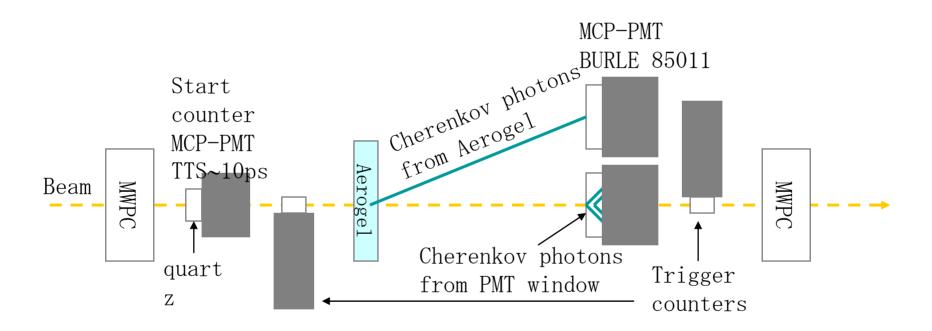



K/π separation by TOF



- Good performance in lower momentum region
- Enable PID under threshold P_c of aerogel

The 3rd SuperB workshop, SLAC

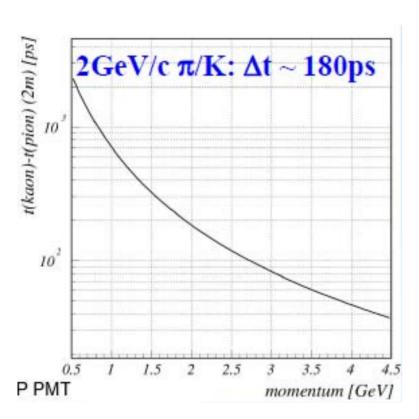

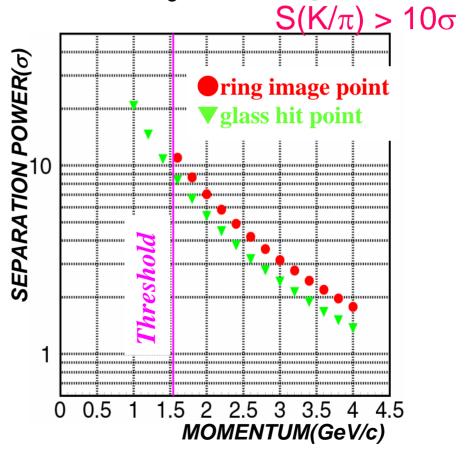

TOP counter MC





Expected performance

bi-alkali photocathode:

 π/K separation at 4GeV/c < 4σ

(because of chromatic

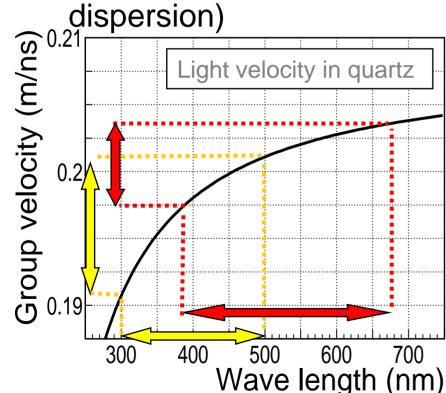

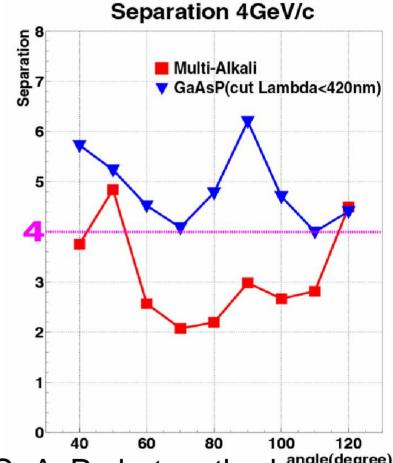


GaAsP photocathode: degree)

 π/K separation at 4GeV/c > 4s less dispersion

TTS



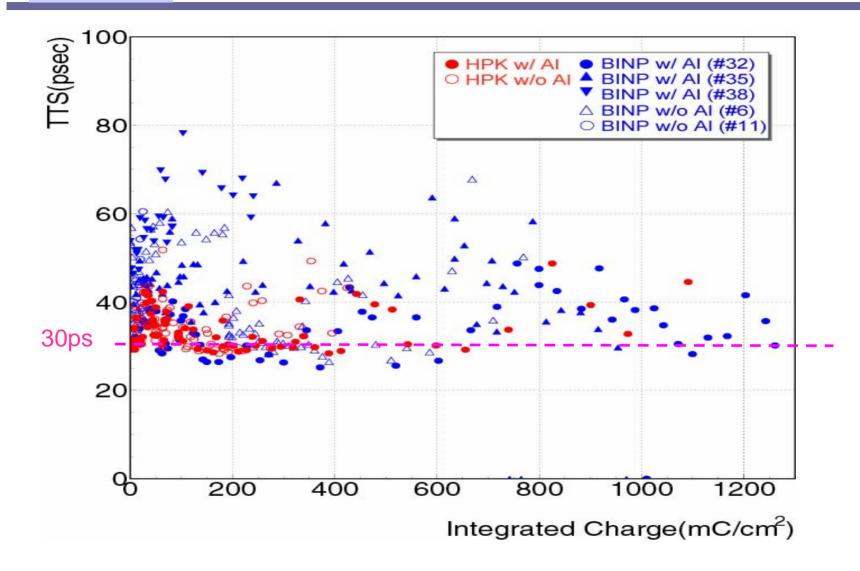

Cross-talk of MCP-PMT




SL10: cross-talk problem solved by segmenting electrodes at the MCP



R&D of Readout ASIC for TOP



- Time-to-Analog Converter → Time resolution <~20ps
- Double overlap gates \rightarrow Less dead time (~100ns)
- 0.35µ CMOS process
- 2nd batch TAC-IC was submitted to VDEC (U. Tokyo)

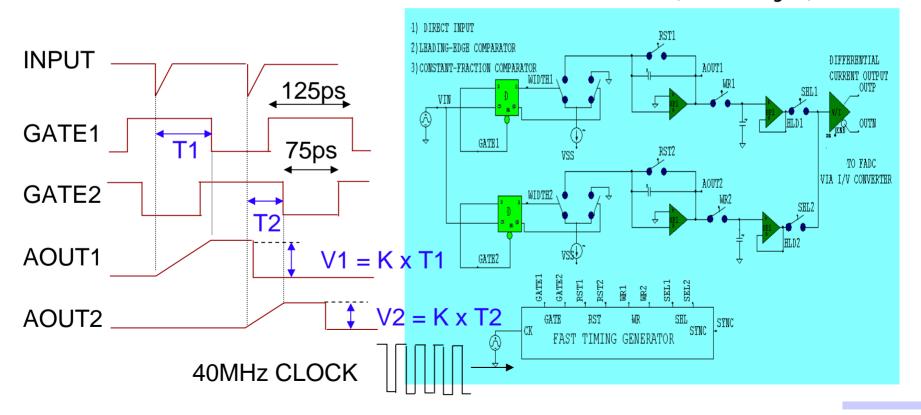

The 3rd SuperB workshop, SLAC

Readout Electronics




Aerogel RICH readout

- Total ~ 100k channels!
- Readout scheme → pipeline
 - Only record hit information


Basic parameters for the ASIC

- **CMOS-FET**
- Gain=10V/pc
- Shaping time=0.15µsec
- VGA=1.25~20
- 18 channels/chip
- Power consumption: 5mW/channel

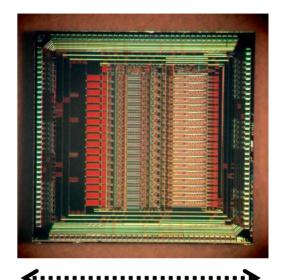


 \square 4.93[mm]

3rd batch was submitted to VDEC (More protection to noise was done)

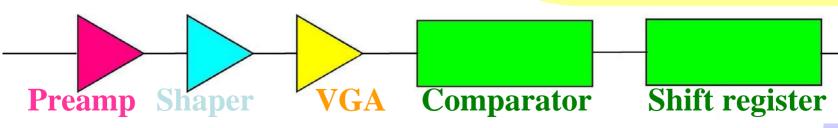


Design


- Quartz: 255cm^L x 40cm^W x 2cm^T
 - cut at θ =46deg. to reduce chromatic error
- Multi-anode MCP-PMT
 - Linear array (5mm pitch), Good time resolution (<~40ps)
 - Three readout plane

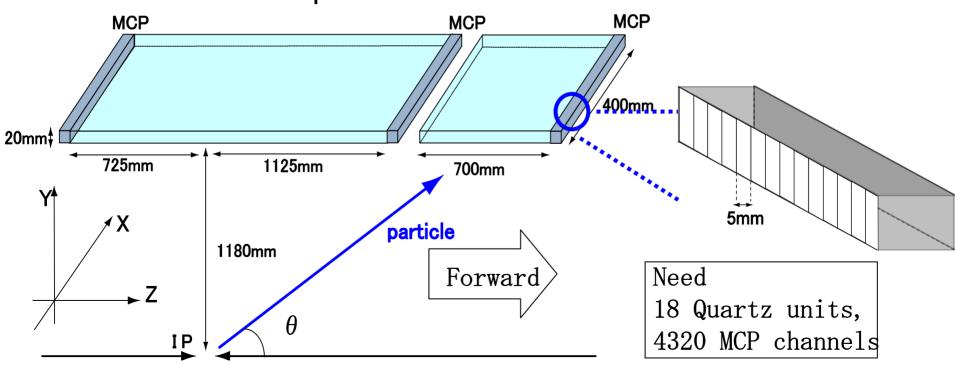


Mechanical design



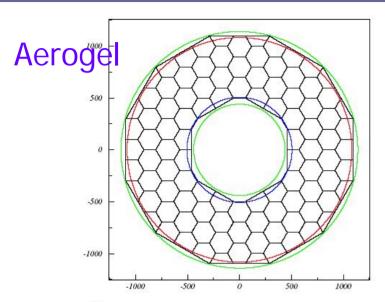

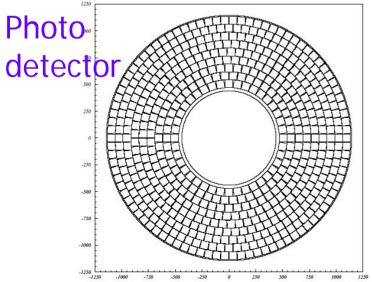

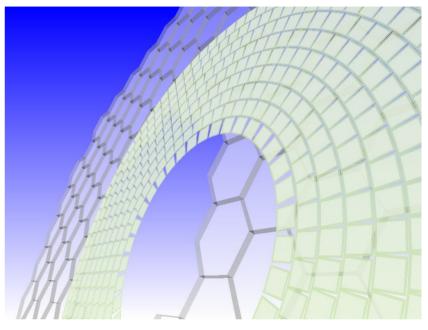

Aerogel radiator

- Hexagonal tiling to minimize aerogel boundary
- side length, 125 mm

Photo detector

- Total PD : 564, 6 sectors
- Cover 89.0% of area

Collaborator



I.Adachi, K.Fujita, A.Gorisek, T.Fukushima, D.Hayashi, T.Iijima, K.Inami, T.Ishikawa, H.Kawai, Y.Kozakai, P.Krizan, A.Kuratani, T.Nakagawa, S.Nishida, S.Ogawa, T.Ohshima, R.Pestotnik, T.Seki, T.Sumiyoshi, M.Tabata, Y.Unno